

ProDetect® Multi-Drug X (2-18) Drugs Rapid Test Cup with/without Adulterant

🎯 High sensitivity, specificity and accuracy. ✅ MDA registered.



Manufactured by:



Medical Innovation Ventures Sdn. Bhd. (MEDIVEN®)
1st Floor, Plot 88F, Lintang Bayan Lepas 10, Bayan Lepas Industrial Park, Phase 4,
11900 Bayan Lepas, Pulau Pinang, Malaysia.
Email: sales@mediven.com.co | export@mediven.com.co
Office / Fax: +604 305 2730
Website: www.mediven.com.my



ABOUT OUR TEST KIT

A visual one-step lateral flow chromatographic immunoassay for the qualitative detection of multiple substances in human urine. Unlike other laboratory diagnostic methods that require complex procedure, equipment and long turnaround time, the test detects drug of abuse substance with a simple test procedure in 5 minutes. The test utilizes monoclonal antibodies to selectively detect elevated levels of specific drugs in urine. The test provides qualitative result visually based on the formation of colored lines, thereby allowing easy interpretation of the test result.



FAST
Test result in 5 minutes.



COST EFFECTIVE
Cost effective method for drug screening.




RELIABLE
High sensitivity, specificity, and accuracy.



TAMPER PROOF
With adulterant test to prevent sample manipulation.

PERFORMANCE CHARACTERISTICS

Diagnostic sensitivity and specificity

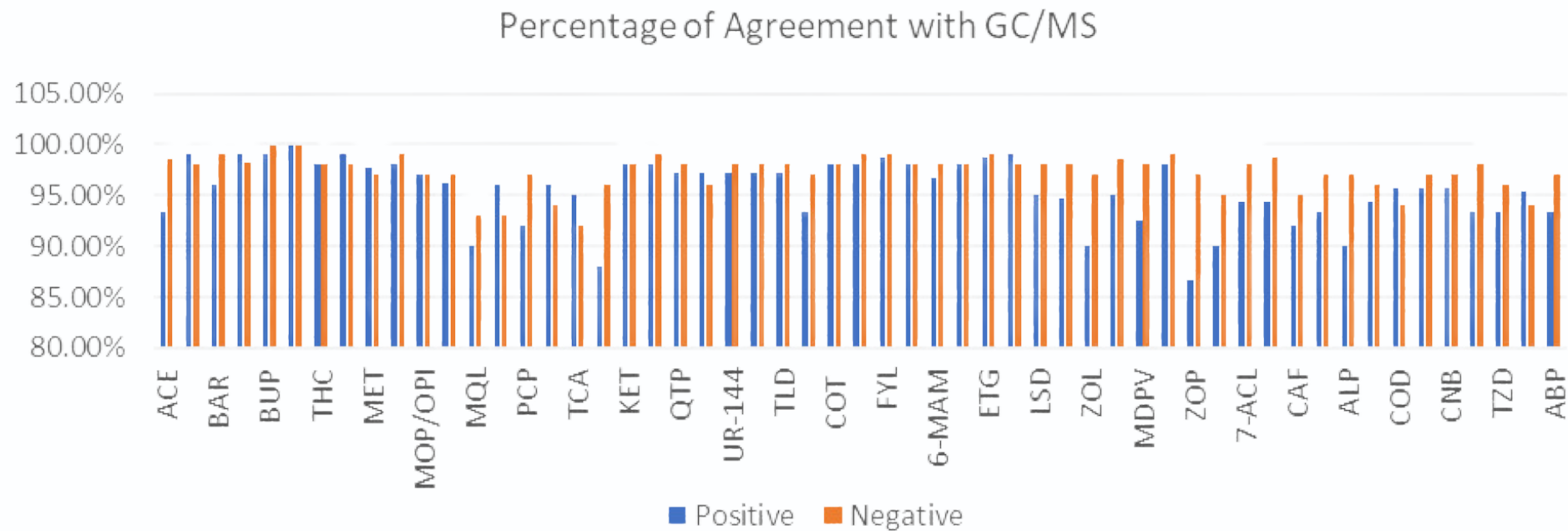

ORDERING INFORMATION

Catalog No.	Product Name	No. of Test Cup/Kit
PR-MDX	ProDetect® Multi-Drug X (2-18) Drugs Rapid Test Cup with/without Adulterant	25

TESTS AND CONCENTRATION

CUT-OFF AVAILABLE

Test	Cut-off (ng/mL)	Test	Cut-off (ng/mL)	Test	Cut-off (ng/mL)	Test	Cut-off (ng/mL)
Acetaminophen (ACE)	5000	Methamphetamine (MET)	1000	Ketamine (KET)	1000	Lysergic Acid Diethylamide (LSD)	10
Amphetamine (AMP)	1000		500		500		20
	500		300		300		50
Barbiturates (BAR)	300	Methylenedioxymethampheta mine (MDMA)	300	Oxycodone (OXY)	100	Methylphenidate (MPD)	300
	200		500		300		150
			1000		100	Zolpidem (ZOL)	50
Benzodiazepines (BZO)	500	Morphine/Opiate (MOP/OPI 300)	300	Cotinine(COT)	300	Mephedrone(MEP)	500
	300		200		200		100
	200		100		100	3, 4- methylenedioxypyrovaleron e (MDPV)	1000
	100	Methaqualone (MQL)	300	2-Ethylidene-1,5- dimethyl-3,3- dipheylpyrrolidine (EDDP)	300		500
Buprenorphine (BUP)	10	Meperidine (MPRD)	100	Fentanyl (FYL)	100	Diazepam (DIA)	300
	5	Opiate (OPI)	2000		300		200
Cocaine (COC)	300		1000		100	Zopiclone (ZOP)	50
	200	Phencyclidine (PCP)	50		20	Methcathinone (MCAT)	500
	150		25		10		300
	100	Propoxyphene (PPX)	300	Synthetic Marijuana (K2)	50	7-Aminoclonazepam (7-ACL)	200
Marijuana (THC)	300	Tricyclic Antidepressants (TCA)	1000		30		100
	200		500		25	Prégabaline (PGB)	50000
	150		300	6-Monoacetylmorphine (6-MAM)	10		500
	50	Tamadol (TML)	100	(+) 3,4-Methylenedioxy- Amphetamine (MDA)	500	Codeine (COD)	200
	30		200	Ethyl-β-D-Glucuronide (ETG)	1000	Gabapentin (GAB)	2000
	25		300		500	Zaleplon (ZAL)	100
	20		500		300	Carisoprodol (CAR)	2000
Methadone (MTD)	300	Carfentanyl (CFYL)	500	Clonazepam (CLO)	400	Alpha- Pyrrolidinovalerophenone (α -PVP)	2000
	200	Cannabinol (CNB)	500		150		1000
Caffeine (CAF)	1000	Tropicamide (TRO)	350	Quetiazepam (QTP)	1000		500
Cathine (CAT)	150	Alprazolam (ALP)	200	Fluoxetine (FLX)	500		300
AB-PINACA (ABP)	10	UR-144	25	Kratom (KRA)	300	Tilidine (TLD)	50